


Universal oven

Model UN450plus

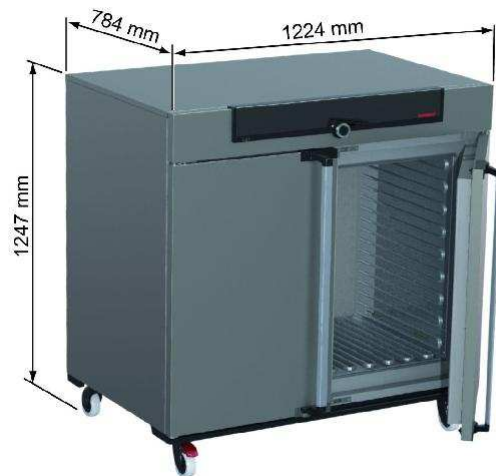

Standard equipment

Ventilation and Control

- natural convection
- adjustment of pre-heated fresh air admixture by air flap control in 10% steps for each segment individually
- vent connection with restrictor flap
- adaptive multifunctional digital PID-microprocessor controller with 2 high-definition TFT-colour displays
- self-diagnostics for fault analysis
- 2 Pt100 sensors Class A in 4-wire-circuit, mutually monitoring and taking over the performance at the same temperature value
- digital timer adjustable from 1 minute to 99 days, 23 hours
- parameters adjustable on ControlCOCKPIT: temperature (Celsius or Fahrenheit), air flap position, programme time, time zones, summertime/wintertime
- SetpointWAIT function - the process time does not start until the set temperature is reached
- HeatBALANCE function for adapting the distribution of the heating performance of the upper and lower heating circuit from -50% to +50%
- Setting of languages on ControlCOCKPIT: German/English/Spanish/French
- resolution of display for setpoint values 0.1 °C up to 99.9 °C, 0.5 °C from 100 °C and for actual values 0.1 °C (LED)
- calibration facility directly at ControlCOCKPIT for three freely selectable temperature values
- programme stored in case of power failure

- incl. works calibration certificate for +160 °C

Multiple Overtemperature Protection

- with audible and visual alarm
- overtemperature monitor TWW, protection class 3.1 or adjustable temperature limiter TWB, protection class 2, selectable on display
- mechanical temperature limiter TB, protection class 1 according to DIN 12 880
- additionally integrated over- and undertemperature monitor "ASF", automatically following the setpoint value at a preset tolerance range; alarm in case of over- or undertemperature, heating is switched off in case of overtemperature
- mechanical temperature limiter TB, protection class 1 according to DIN 12 880 to switch off the heating approx. 10 °C above nominal temperature

Textured Stainless Steel Casing

- w x h x d: 1224 x 1247 x 784 mm
- fully insulated stainless steel doors with 2-point locking (compression door lock)
- rear zinc-plated steel
- on lockable castors

Interior - Heating Concept

- w x h x d: 1040 x 720 x 600 mm, 449 l
- easy-to-clean interior, made of stainless steel, reinforced by deep drawn ribbing with integrated and protected large-area heating on four sides
- 2 stainless steel grids

Temperature Range

- min. 5°C above ambient up to +300°C

Voltage / Power Rating

- 400 V (+/- 10%), 3 phases/N, 50/60 Hz
- ca. 5.800 W

Packing Data

- net weight approx. 170 kg
- gross weight carton approx. 240 kg
- dimensions approx.:
carton w x h x d: 133 x 145 x 105 cm
- the appliances must be transported upright

Customs Tariff Number

- 8419 8998

Country of Origin

- Federal Republic of Germany

WEEE-Reg.-No.

- DE 66812464