




circulation baths

WWW.FALCINSTRUMENTS.IT

circulation baths

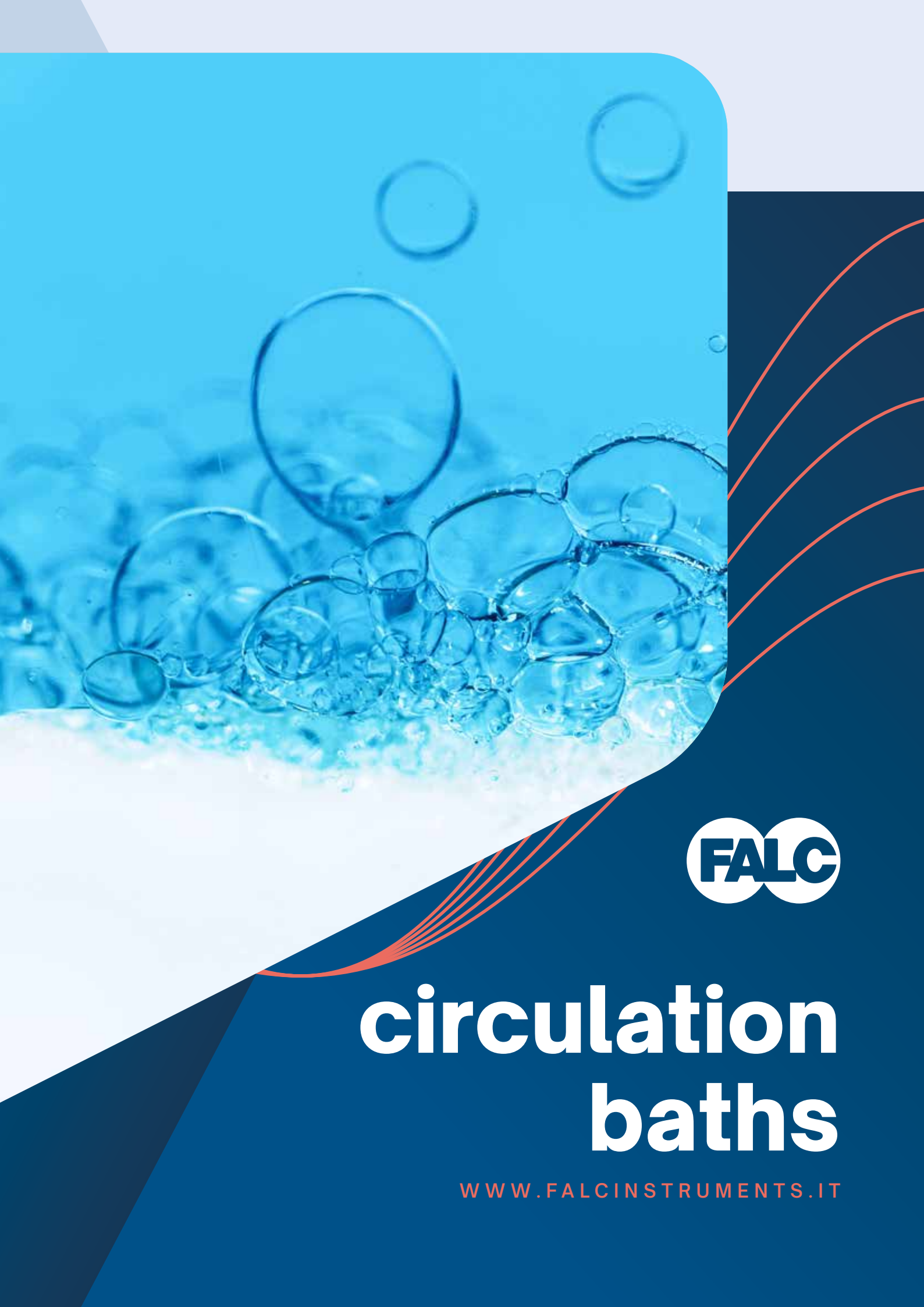



I N D E X

04

circulation baths

Circulation-heating baths

05

Circulation-cooling baths

circulation baths

CIRCULATION-HEATING BATHS



Scan for
more info
about the
product



SB



LCD touchscreen display
The tank is made of 18/8 stainless steel pressed in one piece with lid



The temperature is set and checked by a digital thermostat with PID controller and PT100 probe. The stirring pump ensures excellent temperature uniformity inside the tank
Digital timer hh:mm:ss or infinite



The instrument is designed for connection to the external pump set
Manual operation of the fireplace
Setting for 4 languages (Italian-English-French-Spanish)
Program saving function (temperature/maintenance/pump speed)
Programmable delayed start (day hour)

Technical data		SB 5T	SB 15T	SB 24T
Capacity tanks	lt	5	12	18,5
Drain tap		no	no	si
Digital timer	hh:mm:ss	99:59:59 or ∞	99:59:59 or ∞	99:59:59 or ∞
Adjustable temperature	°C	+5 amb. to 150	+5 amb. to 150	+5 amb. to 150
Precision with H ₂ O	°C	±0,5	±0,5	±0,5
Display resolution	°C	0.1	0.1	0.1
Pump flow rate mbar 190	lt/min	9	9	9
Head m. 1.8	lt/min	7,8	7,8	7,8
Heating element power	W	1300	1300	1300
Internal dimensions LxWxH	mm	250x100x150	290x255x150	240x440x150
Usable tank dimensions LxWxH	mm	100x125x150	255x165x150	255x330x150
External dimensions LxWxH	mm	330x180x450	350x360x450	350x530x450
Weight	kg	7	8	11
Code		612.1029.06	612.1029.16	612.1029.25

V/HZ 220-240/50-60

SAFETY CLASS 0

PROTECTION CLASS IP 30


ACCESSORIES



EXTERNAL PUMP SET
It enables the external power supply via a closed circuit.
Ø12mm tube.
Code 612.1029.02



Stainless steel floating level device for water shortage
Code 612.1029.04

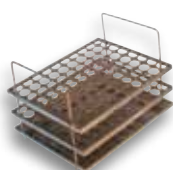

18/8 stainless steel tube rack with 3 shelves

	1 per tank	1 per tank	2 per tank
	12 holes Ø21	48 holes Ø21	48 holes Ø21
Code	212.1250.02	212.1250.05	212.1250.05
	20 holes Ø17	60 holes Ø17	60 holes Ø17
Code	212.1250.03	212.1250.06	212.1250.06

circulation baths

CIRCULATION-COOLING BATHS



LCD touch screen display
Tanks are made of one piece
of 18/8 stainless steel sheete



The temperature is checked by a PID microprocessor
and PT100 probe



Automatic shutdown of the compressor over 37.0°C



Digital timer hh:mm:ss or infinite



Eco Friendly function, the cooling compressor starts working
only when necessar



Easy access to keep the instrument clean and effective at all times
Manual operation of the fireplace
Setting for 4 languages (Italian-English-French-Spanish)
Program saving function (temperature/maintenance/pump speed)
Programmable delayed start (day hour)



SBF

Technical data		SBF 7T	SBF 15T
Capacity	lt	7	15
Drain tap		no	no
Timer	hh:mm:ss	∞ or 99:59:59	∞ or 99:59:59
Working temperature	°C	-30 to +150	-30 to +150
Precision in the tank	°C	±0.2	±0.2
Heating capacity	W	1300	1300
Cooling capacity	°C	-25 0 25	-25 0 25
	W	210 170 260	520 600 850
Pump type		Pressure	Pressure
Pump Flow rate	lt/min	7.8	4
Prevalence	mt	2.5	2.5
Bottom dimensions LxWxH	mm	100x250x200	250x290x200
Usable tank size LxWxH	mm	100x150x200	250x190x200
Outer dimensions LxWxH	mm	300x440x730	440x475x730
Weight	kg	31	41
Code		612.1050.08	612.1050.16

V/HZ 220-240/50-60

SAFETY CLASS 2

PROTECTION CLASS IP 30

ACCESSORIES



Stainless steel floating level device for water
shortage

Code 612.1029.04

Wheels

Code 612.1050.31

External circulation pipes 1.5 m

Code 612.1050.33

1 lt liquid for cryostat -30÷80°C 1 LT

Code 612.1050.35

Notes

[illegible]

Notes

[illegible]



FALC INSTRUMENTS s.r.l.

Via G. M. Compagnoni, 2
24047 Treviglio (BG) - Italy
+39 0363 304660 | falc@falcinstruments.it

WWW.FALCINSTRUMENTS.IT