


# CONTROL PANELS

## LGK-2 Digital

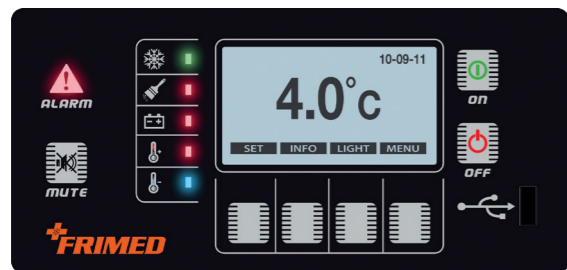

Microprocessor control panel with easy-to-read **digital display** showing inside menu, temperature, set point, date and hour, model of appliance.

Equipment settings and switching off only by means of **password**.

Continuous **temperature monitoring** with set point and alarms shown on display.

**Visual and acoustic alarms** through buzzer sound, general alarm warning light, specific alarm warning light and icon on display.

Alarms included: high and low temperature (adjustable threshold), open door, condenser efficiency, damaged probes, power failure, run-down battery.

**Back up battery** (optional) to show temperature on display and generate alarms in case of power failure.

**Remote contact** to connect the appliance to an external alarm.

**GSM module** (optional) to send SMS alarm messages showing date, time, cabinet temperature and type of alarm. Possibility to send SMS message to check correct functioning of cabinet.

**LAN port** (optional) to connect one or more appliances via ethernet on PC for monitoring correct functioning.

**USB port** to download temperature data and alarms or to reset the appliance.

**Data logger** to record temperatures of last 30 days. Download data on USB pen drive.

**Memory** for last 20 alarms shown on display.

