

Key switch

Digital temperature display

USB-Port for reading out data

Minimum/maximum temperature memory

Visual and audible alarm signal

Comfortable access


potential-free contact/  
RJ 45 socket by optional  
PC-KIT-NET

Antifreeze

Forced-air cooling

5 Shelves

Automatic fast defrosting with removable condensation container

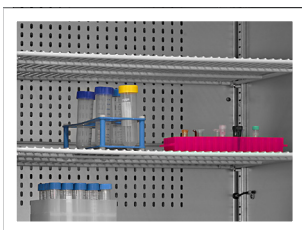

**Interior protected**

Possible static charges are prevented by a dissipation system.



**Explosion protection - Zone 2**

The interior is free from ignition sources.

# LABEX® 520/720

## PRO-ACTIVE



### PRO-ACTIVE- Control:

Permanent, proactive monitoring of the performance data and alerts in the case of deviations – so that you can take countermeasures in good time before a fault threatens your chilled goods; world's most accurate temperature control in refrigerated areas thanks to two standard PT-1000 sensors.

- **External housing** made from galvanised sheet steel (rust-proof), with white, anti-scratch powder coating. Length of the plug cable: approx. 3 m.
- **Equipped with castors.**
- **Interior** made from smooth aluminium with colourless protective coating. Gibs (every 15 mm) for flexible interior arrangement.
- **Interior consists** of five shelves on brackets.
- **Extra-thick 70-mm insulation**, made from high-quality, compression-moulded and environmentally friendly material. Energy saving.
- **Self-closing door** with easy-to-replace plastic magnetic seal frame, lockable.
- **Door hinge** on the right-hand side as standard (see illustration), or the left-hand side, can be retrofitted.
- **Forced-air cooling**, switches off automatically when you open the door, ensures a uniform temperature and minimizes temperature deviations.

- **Valve controlled pressure compensation** for easier door opening.
- **Automatic fast defrosting** with time limit and temperature monitoring. Defrost sensors are protected by safety barriers.
- **Removable condensation container.**
- **Power failure alarm** (visual and audible alarm), the monitoring unit remains in operation for approx. 30 hours on battery power.
- **Key switch** protected control panel against tampering. Membrane keyboard with digital temperature display and minimum/maximum temperature memory. Control and display sensors protected by safety barriers.
- **Warning functions** with visual and audible alarm signal in the case of temperature deviations or other malfunctions. Door open alarm after 60 seconds.
- **Alarms** can be forwarded using a potential-free contact (e.g., to a mobile phone with optional GSM-Module or to a control centre).
- **Data documentation** can be read via USB interface with KIRSCH Datanet software.
- **Antifreeze** against sub-zero temperatures.
- **Ventilation-enforced refrigerating machine**, hermetically sealed, energy saving, low-noise, easy to service.

## Specification

	LABEX® 520 PRO-ACTIVE	LABEX® 720 PRO-ACTIVE
Capacity	500 Liter	700 Liter
Temperature setting	approx. +0 °C to +15 °C	approx. +0 °C to +15 °C
Voltage	220 – 240 V, 50 Hz	220 – 240 V, 50 Hz
Power consumption	234 Watt	237 Watt
Normal consumption	1,2 kWh/24 h	1,23 kWh/24 h
Admissible ambient temperature	from +10 °C to +38 °C	from +10 °C to +38 °C
Heat emission (max.)	543 Watt	545 Watt
Exterior dimensions (including door handle and distance from wall)	w x d x h = 77 x 76 x 193-196 cm	w x d x h = 77 x 98 x 193-196 cm
Interior dimensions	w x d x h = * 62 x 57 x 140 cm	w x d x h = * 62 x 77 x 140 cm
Exterior dimensions with door open at 90°	w x d = 77 x 144 cm	w x d = 77 x 165 cm
Shelfsize	w x d = 59 x 45 cm	w x d = 59 x 65 cm
Max. load shelf	40 kg	40 kg
Weight	netto 120 kg, brutto 148 kg	netto 145 kg, brutto 174 kg

\*usable width 2 cm, usable depth 11 cm, usable height 13 cm less.

## Optional equipment

- **Glass door with lock.**
- **Adjustable feet** to compensate uneven floors.
- **Drawer conductive on roller runners**  
Only for LABEX® 520 PRO-ACTIVE: w x d x h = 57 x 39 x 5.6 cm (a maximum of twelve drawers is possible).
- **Additional length and cross dividers** (for drawer).
- **Additional shelf** on rails or brackets.
- **Basket on rails**  
LABEX® 520 PRO-ACTIVE: w x d x h = 60 x 45 x 10 cm, LABEX® 720 PRO-ACTIVE: w x d x h = 60 x 65 x 10 cm.
- **Aluminium tray with rails**  
LABEX® 520 PRO-ACTIVE: w x d x h = 60 x 40 x 2.5 cm, LABEX® 720 PRO-ACTIVE: w x d x h = 60 x 65 x 2.5 cm.
- **Outer housing in stainless steel 4301**, longitudinal brushed.
- **GSM-Module** for alarm text message transmission (e.g., to a mobile phone).
- **Water-cooled refrigerating machine.**
- **Temperature documentation.**
  - Temperature sensor PT 100 or PT 1000, **with safety barriers** and reference body (measuring bottle or cold block) is installed.
  - 4-pin PT 100 or PT 1000 temperature sensor, class 1/3 B, including built-in reference body.
  - Mechanical pen-recording thermometer with two waxed paper strips.