


orto alresa

Centrifuges and laboratory products

info@ortoalresa.com

www.ortoalresa.com

Tel.: +34 91 884 40 16

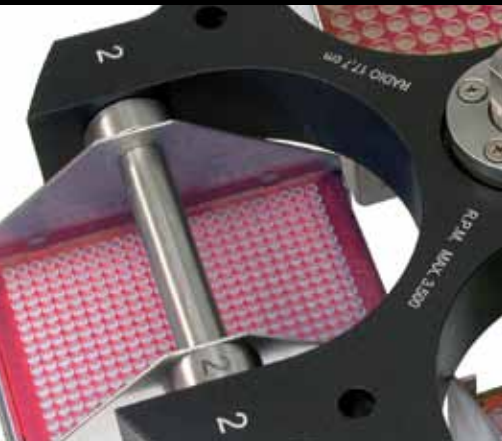

Special applications: Centrifuges: Life sciences

VETCEN

VETCEN

Small veterinaries need versatile equipment, friendly to use and efficient. The centrifuge **Vetcen** has two rotors which allow run samples on Capillaries tubes, microtubes of 1,5 to 2,2 and 5, 10 and 15 ml tubes.

The versatility of this equipment is given by its multi rotor tubes of different volume can in one run , spin capillary tubes and microtubes. Optimal for processing a small number of samples is undoubtedly the most appropriate choice for emergency laboratories, and getting shortened preparation time, simplifying the separation of samples.

User friendliness

- LED screen shown (RPM on steps of 100, RCF on steps of 10xg) time and deceleration.
- Intuitive use with start, stop, lid opener and short centrifugation buttons. Blocking/modifying option of RPM/RCF while operating.



- Timer count up/down, from 0 or at "set RPM/RCF" for test reproducibility. Rotor list on memory.
- Last values remain on memory. Acceleration and braking selectable by operator.
- Microprocessor controlled. Induction motor maintenance free.
- Automatic lid latch lock. Open lid automatically selectable.

Features

- Capacity: Max. 6 x 1.5-2.2 ml + 6 x 1.5 x 75 ml.
- Speed max. 11.500 RPM/12.716 xg.
- Noiseless <60 dB.

Safety

- Lid locking and holding.
- Manual overriding. Lid dropping protection. Imbalance switch off.
- Chamber of centrifugation in stainless steel. Adapters and rotors autoclavables.
- Viewing and calibration port in the lid. Guard ring.
- Several messages keep the operator informed about the device situation.

EU Directives: 98/79/EC, 2006/95/EC, 2004/108/EC, 2011/65/EU, 2012/19/EU.
Standards: EN 61010-1, EN 61010-2-020, EN 61326-1.

Versions

Code	Dimensions (mm)(w x d x h)			Net weight (Kg)	Voltage (V)	Frecuency (Hz)	Consumption (W)
CE 160	276	390	272	16	220-240	50-60	140
CE 177	276	390	272	16	110-120	50-60	140

Accessories



ROTOR	ANGLE FIXED 30 °		MIXTO		
	Max. capacity	12x1,5x75 mm.	6x1,5x75+6x1,5/2 ml.		
RPM	11.500		11.500		
Radius(mm)	86		86		
RCF Max. (xg)	12.716		12.716		
SAMPLE VOLUME	Dim (mm) approx.	ADAPTERS		ADAPTERS	
		Tubes	Ref.	Tubes	Ref.
1.5-2.2 ml	Ø 11x39	-	-	6	-
Capillaries	1,5 x 75 mm.	12	-	6	-

orto alresa

About centrifugation

info@ortoalresa.com

www.ortoalresa.com

Tel.: +34 91 884 40 16

Álvarez Redondo S.A.
C/ Misericordia, 23.
Fax: + 34 91 884 44 23
28864 AJALVIR Madrid (Spain)
Edition: january 2014

