

The professional electric pipetting solution

WITOPED XP is a professional electric pipetting aid, suitable for all pipettes up to 100 ml. The special features like speed adjustable suction and dispensing buttons, the comfortable base which allows stable placing of the device inclusive attached pipette as well as the high-precision, long lasting motor make **WITOPED XP** an indispensable device in your laboratory.

Nosepiece release

This button protects against accidentally releasing the nosepiece.

Removable stand

The removable stand allows stable placing of the device along with attached pipette.

Nosepiece

The nosepiece is suitable for all kinds of pipettes. In case of contamination, the nosepiece can be easily removed in order to exchange the filter and the silicone adapter.




Dispenser buttons

The buttons can be individually adjusted to your application by turning:

- (F) Fast = 10 ml/s
- (M) Medium = 5 ml/s
- (S) Slow = 1 ml/s
- (G) Gravity drain

Motor

The high-precision motor is robust and made for a long life time.

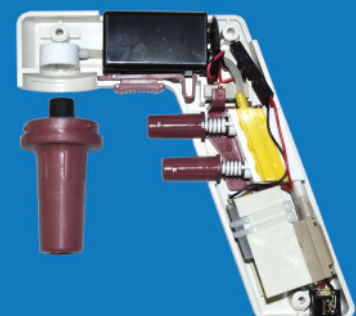

Pipetting aid

Your advantage at a glance:



- high-quality instrument, most competitive price
- high-precision, durable motor
- individual adjustment of suction and dispensing speed
- three pipetting speeds selectable:
 - (F) Fast = 10 ml/s
 - (M) Medium = 5 ml/s
 - (S) Slow = 1 ml/s
- (G) Gravity drain
- gravity drain mode for draining out without using any motor
- handy and light weighted for fatigue-proof pipetting
- tight fit of any pipette into the pipette socket
- nosepiece with replaceable filter and silicone adapter
- nosepiece autoclavable at 121°C
- removable stand for stable holding of instrument along with attached pipette
- including wall bracket
- including NiMH 8,4V 200mAh battery
- CE label

WITOPED XP Professional

Type	Order number
220 V	5 380 001
120 V	5 380 002


Scope of supply

WITOPED XP, battery charger, replacement filters, removable stand, instruction manual



Replacement one-way filter

Packaging unit	Order number
5	5 380 006

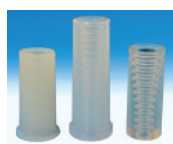

Replacement silicone adapter

Order number
5 380 004