


“witeg” **WiseStir® HT-AX** High-Torque Analog Overhead Stirrer, Max.40-/60-Lit., Max.1.000rpm with Chrome-Plated & Stainless Steel Chuck , Easy-to-Use Dial for Speed Setting, for High-Viscosity Fluids, 150.000 mPas

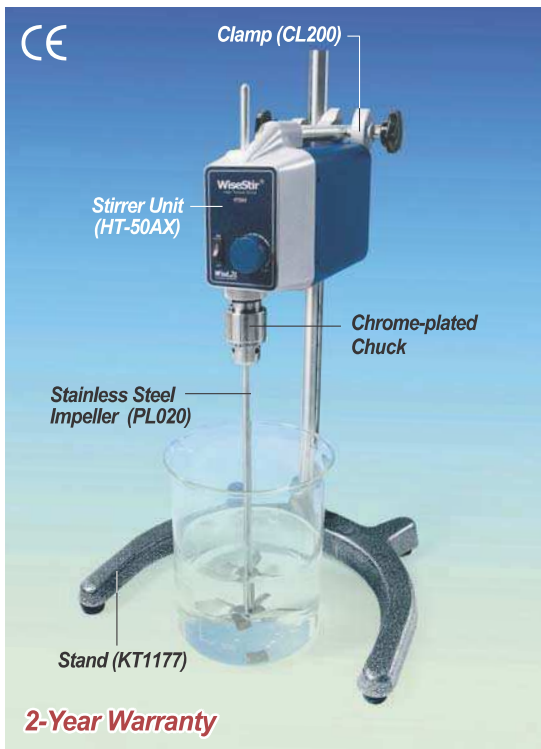

2-Year Warranty

High-Torque Analog Overhead Stirrer,
□HT-50AX□Package-set

- High Torque Stirring by 3,5:1 Ratio Worm-gear Motor
- Smooth Start and Change of the Stirring Speed
- CE Certified

- * Ideal for Stirring Material of Viscosities Ranging from Light Watery Solutions to High-Viscosity Fluids
- * Ideal for Schools, Universities and Inspection Laboratories
- * Food, Beverage, Dairy, Sugar & Candy Production
- * Smooth Start and Change of the Stirring Speed
- * Easy-to-use Dial for Speed Setting
- * The Position(height) of Impeller is Adjustable
- * Over-Heat Protection Function for Motor
- * Chrome-Plated & Stainless Steel Chuck-Prevention from Corrosion

Specification

Stirring Capacity & Models	Max. 40 Lit. (HzO) □HT-50AX□	Max. 60 Lit. (HzO) □HT-120AX□
Maximum Viscosity	100.000 mPas	150.000 mPas
Maximum Torque	300 Ncm	500 Ncm
Speed Range & Motor	0 ~ 1.000 rpm, Worm-Geared DC Motor (3,5:1 Ratio)	
Motor Rating (Input / Output)	83 / 75 W, 1 / 10 HP	130 / 120 W, 1 / 6 HP
Controller	Analog Phase Controller	
Max. Impeller Shaft Diameter	10 mm	
Permissible Temp. & Humidity	0 ~ 40°C, 80%	
Safety Device	Over Heat Protector for motor	
Material	ABS and Aluminum	
Dimension & Weight	90×180×h130mm, 3,5 kg	90×210×h130mm, 4 kg
Packing Size & Gross Weight	320×300×h180mm, 5 kg	320×300×h180mm, 8 kg
Power Consumption	85 W	130 W
Power plug type	Hybrid E+F	
Power Supply*	AC 230V, 1N~, 50/60Hz	

* Other Specifications are available upon Customer s Request.

Ordering Information

- * □HT-50AX□Package-set : Complete with ① Overhead Stirrer Unit ② □PL020□Stainless Steel Impeller (Rod φ8×L500mm, Blade φ90mm) ③ □CL200□Clamp ④ □KA1177□Safety Stand (Rod φ23×L800mm, Stand Base w34×d29cm)
- * □HT-120AX□Package-set : Complete with ① Overhead Stirrer Unit ② □PL020□Stainless Steel Impeller (Rod φ8×L500mm, Blade φ90mm) ③ □CL200□Clamp ④ □KA1177□Safety Stand (Rod φ23×L800mm, Stand Base w34×d29cm)

With 40 Lit. Stirring Capa., 100.000 mPas, “HT-50AX” Package-set & Unit

DH.WOS01252 High Torque-type Overhead Stirrer, **Package-set, 230V**
DH.WOS01251 , **Unit only, 230V**

With 60 Lit. Stirring Capa., 150.000 mPas, “HT-120AX” Package-set & Unit

DH.WOS01272 High Torque-type Overhead Stirrer, **Package-set, 230V**
DH.WOS01271 , **Unit only, 230V**

※ You will find matching beakers in our general catalogue at following pages : 292

※ Refer to page 156, 157 for Accessory Specifications.

※ Refer to page 149 for Product Line & Selection Guide



Analog High-Torque Overhead Stirrer,
□HT-120AX□Unit only